




电话: 0755-84870203

邮箱: sales@highlightoptics.com





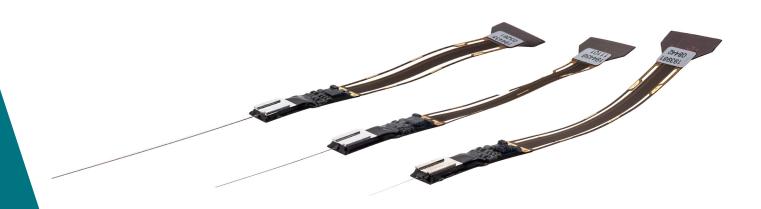

NEUROPIXELS 1.0 NHP SHORT-MEDIUM-LONG

High resolution fully integrated silicon neural probe for large animal recording

Key Features

- 960, 2496 or 4416 reliable, low-impedance TiN electrodes
- Dense 2 row linear electrode layout along one 10 mm, 25 mm or 45 mm long single shank
- 125 µm wide (NHP MEDIUM and NHP LONG) or 70 µm wide (NHP SHORT) and ~110 µm +/- 15 µm thick shank
- Maximal shank bending ≤200 µm
- 384 parallel, dual-band (AP¹, LFP²), low-noise recording channels
- On-chip amplification, signal conditioning and digitization

- Channel-independent configuration and reference selection (internal or external)
- Small, flexible and light-weight package (0.4 g)
- Systematic quality control process to ensure low variability in performance
- Compatible with SpikeGLX and Open Ephys software
- Compatible with the Neuropixels 1.0 headstage, cabling and PXIe system or OneBox


Important Information

The Neuropixels probes are intended for RESEARCH USE ONLY ("RUO") in non-human subjects such as small animals*. These Neuropixels probes should not be used in humans and are not manufactured or approved for human use. They have no proven human efficacy and are not indicated for human use or any form of clinical use. The Neuropixels probes are provided and delivered for use only under the imec general terms and conditions of sale of Neuropixels 1.0 probes ("GTC"). [The GTC is available for download on www.neuropixels.org]

Action potentials

² Local field potentials

Description

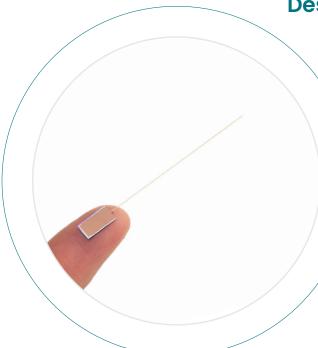

The Neuropixels³ 1.0 NHP neural probe is a CMOS digital integrated microsystems platform, based on the standard Neuropixels 1.0 probe. A thicker and longer shank will allow for high density recording in large animal models.

The probe features 960, 2496 or 4416 low-impedance
TiN recording sites densely tiled along one thin, 10, 25 or
45 mm-long, straight shank. The 384 parallel, configurable,
low-noise recording channels integrated in the base
enable simultaneous, dual band recording of hundreds of
neurons

On-chip circuitry for signal conditioning and digitization results in a small and light-weight package allowing the implantation and simultaneous use of multiple probes in close proximity.

Neuropixels probes enable long-term monitoring and dense sampling of single cell activity as well as larger neuron populations in awake and anaesthetized animals. This probe is compatible with the Neuropixels 1.0 readout system.

3 JJ Jun et al., Nature 2017, 551, 232-236

Ordering information



ORDER CODE	DESCRIPTION
NP1015	Box of 5 Neuropixels 1.0 NHP short 10mm probes
DNP1015	Box of 6 Neuropixels 1.0 NHP short 10mm dummy probes
NP1022	Box of 5 Neuropixels 1.0 NHP medium 25mm probes
DNP1022	Box of 6 Neuropixels 1.0 NHP medium 25mm dummy probes
NP1032	Box of 5 Neuropixels 1.0 NHP long 45mm probes
DNP1032	Box of 6 Neuropixels 1.0 NHP long 45mm dummy probes
HS_1000	Headstage for Neuropixels 1.0 probes
HOLDER_1000_C	Neuropixels 1.0 metal cap probe holder pair

Key Applications

- High-density in vivo recording of neural activity in larger animal models such as non-human primates.
- Recording of large neuron populations from several brain regions in freely moving animals at high spatiotemporal resolution and large volume coverage.

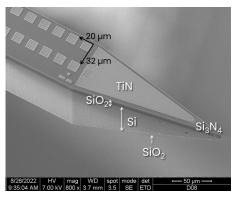


Figure 1: SEM image of the NHP SHORT shank tip. Indicated are the electrode pitch and exposed materials.

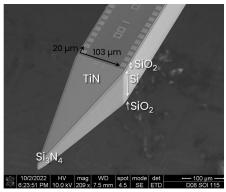


Figure 2: SEM image of the NHP MEDIUM and LONG shank tip. Indicated are the electrode pitch and exposed materials.

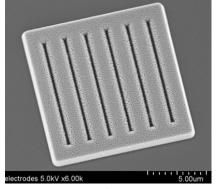


Figure 3: SEM image of a 12 x 12 μ m TiN electrode.



电话: 0755-84870203

邮箱: sales@highlightoptics.com

Key Specifications

ELECTRODES	NHP SHORT	NHP MEDIUM	NHP LONG
NUMBER	960	2496	4416
PATTERN	LINEAR	LINEAR	LINEAR
PITCH	32 µm (column) 103 µm (column) - 20 µm (row) 20 µm (row) (see Figure 1)		
MATERIAL	Porous TiN ⁴ (Figure 3)		
SIZE	12 x 12 µm		
IMPEDANCE	~150 kΩ (at 1 kHz in PBS ⁵)		
SELECTIVITY	Local switch under each electrode		

SH	Λ٨	IK	PR	$\cap P$	FD.	TIFS
OH.	Αľ	<i>1</i> N	Γ	CJE	Γ	LIEO

AND MATERIALS	NHP SHORT	NHP MEDIUM	NHP LONG
NUMBER		1	
WIDTH	70 μm	125	μm
LENGTH	10 mm	25 mm	45 mm
THICKNESS		122 µm	
BENDING	≤100 µm (base to tip)		
TIP LENGTH	175 µm 342 µm		
TIP SHAPE	Sharpened chisel		
TIP ANGLE	~20° in plane and ~25° out of plane (Figure 1 and 2)		
FRONTSIDE MATERIAL	Silicon nitride (Si ₃ N ₄) (Figure 1 and 2)		
BACKSIDE MATERIAL	Silicon dioxide (SiO ₂)		
SIDEWALL MATERIALS	Silicon (Si), silicon dioxide (SiO ₂)		

RECORDING CHANNELS AND

DIGITAL INTERFACE	NHP SHORT	NHP MEDIUM	NHP LONG
NUMBER		384 (dual-band)	
AP BANDWIDTH		0.3-10 kHz	
LFP BANDWIDTH		0.5-500 Hz	
AP INPUT-REFERRED NOISE		5.9 µV _{rms} (typical ⁵)	
LFP INPUT-REFERRED NOISE		9.2 µV _{rms} (typical)	
AP SAMPLING FREQUENCY		30 kHz	
LFP SAMPLING FREQUENCY		2.5 kHz	
DIFFERENTIAL GAINS		50-3000 (8 values)	
CROSSTALK		≤0.13% (at 1 kHz; typical)	
INPUT VOLTAGE RANGE		±5 m _{vpp}	
ADC RESOLUTION		10 bits	
DATA RATE		163.8 Mb/s	
POWER CONSUMPTION		~15 mW (in recording mode; typical)	
SHANK HEATING		<1°C (in the brain)	

REFERENCE SELECTION	NHP SHORT	NHP MEDIUM	NHP LONG	
INPUTS	3 internal	7 internal	12 internal	
	Large tip electrode on the shank (Figure 1 and 2)			
	External input on the probe package (Figure 4)			





电话: 0755-84870203

邮箱: sales@highlightoptics.com

PACKAGE DESCRIPTION	NHP SHORT	NHP MEDIUM	NHP LONG
WIDTH AT PROBE BASE (WI)		6.2 mm	
WIDTH AT SMD ⁷ BASE (W2)		7.2 mm	
WIDTH OF METAL CAP (W3)		4.8 mm	
WIDTH OF FLEX (W4)		4.3 mm	
LENGTH OF PROBE BASE (L1)		10.7 mm	
LENGTH OF SMD ⁷ BASE (L2)		12.2 mm	
LENGTH OF METAL CAP (L3)		7.3 mm	
LENGTH OF FLEX (L4)		39.5 mm	
THICKNESS AT PROBE BASE		~1.8 mm (with metal cap)	
THICKNESS OF FLEX		80 µm	
EXTERNAL REFERENCE INPUT		REF (multiple pads along flex)	
GROUND INPUT		GND (multiple pads along flex)	
BLACK EPOXY		epo-tek / h70e	
CONFORMAL COATING OF SMD ⁷		elpeguard / SL 1307 FLZ-T	
WEIGHT		440 mg (with metal cap)	

HEADSTAGE	NHP SHORT	NHP MEDIUM	NHP LONG
SIZE		15 x 16 mm	
WEIGHT		0.9 g	
ZIF CONNECTOR		45-pin	
CABLE CONNECTOR		4-pin (Omnetics)	
LED INDICATOR		One red LED	
MECHANICAL FIXTURES		Two mounting holes of 1 mm Ø	
CONFORMAL COATING OF SMD ⁷		ELPEGUARD / SL 1307 FLZ-T	

METAL CAP HOLDER	NHP SHORT	NHP MEDIUM	NHP LONG
LENGTH		17.8 mm	
DIAMETER		6.25 mm	
MATERIAL		Aluminium 6061	

⁷ Surface-mount devices: Biasing resistors, decoupling capacitors, EEPROM with probe ID, low-noise reference supply IC

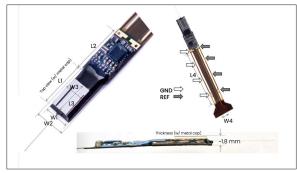


Figure 4: Dimensions of the different probe packages and locations of REF/GND input pads..

The Neuropixels 1.0 neural probe is an advanced silicon CMOS digital integrated microsystem and a tool for neuroscience research. It was developed through a collaboration funded by Howard Hughes Medical Institute (HHMI), Wellcome Trust, Gatsby Charitable Foundation and Allen Institute for Brain Science. Probes were designed, developed and fabricated at imec, Leuven, Belgium in collaboration with HHMI Janelia Research Campus, Allen Institute for Brain Science and University College London.

Legal Disclaimer

The contents of this document are provided by imec, 'as is'. Imec makes no representations nor warranties with respect to the accuracy or completeness of the contents of this publication and reserves the right to make changes to the specifications at any time without notice. All trademarks are the property of their respective owners. The Neuropixels 1.0 probes are ONLY distributed to the Neuroscience Research Community under the "IMEC GENERAL TERMS AND CONDITIONS OF SALE OF NEUROPIXELS 1.0 PROBES". Imec is a registered trademark for the activities of IMEC International (a legal entity set up under Belgian law as a "stichting van openbaar nut"), imec Belgium (IMEC vzw supported by the Flemish Government), imec the Netherlands (Stichting IMEC Nederland, part of Holst Centre which is supported by the Dutch Government), imec China (IMEC Microelectronics (Shanghai) Co. Ltd.), imec India (Imec India Private Limited), and imec Florida (IMEC USA papelectronics design center). nanoelectronics design center).

^{*} Small animals like rodents and non-human primates



电话: 0755-84870203

邮箱: sales@highlightoptics.com